**Tecal 425S**  
**Calibrator**  
Operator's Manual

Issue 4.

Date of issue:

23rd July 2003

® The Techne logo is the trademark of

Techne and Techne Inc.

© Techne , 2003


**Barloworld** Scientific

**CONTENTS****Tecal Calibrator Operator's Manual**

	<i>page</i>
<b>SAFETY AND INSTALLATION</b>	
Declaration of conformity	3
English	4
Français	5
Deutsch	6
Español	7
<b>THE CALIBRATOR</b>	8
Applications	8
Warning	8
Specification	9
Working Conditions	9
Packing	9
<b>WHEN YOU SWITCH ON</b>	10
Front Panel controls	10
Setting the Units for Temperature	12
Setting the Operating Temperature	12
Cooling Probe	13
After Use	13
<b>RS232 SERIAL INTERFACE</b>	14
Sending data to the PC (Auto Control)	14
<b>TECHNICAL INFORMATION</b>	15
General advice	15
General Fault finding	15
Replacement parts	16
Accessories	16
Inserts	17

Techne  
Duxford  
Cambridge  
CB2 4PZ

## Declaration of Conformity

Techne Unit Tecal 425S has been designed to comply with the following European Standards:

EN 50081-1:1992 Electromagnetic Compatibility; Generic emission standard.

EN 50082-1:1992 Electromagnetic Compatibility; Generic immunity standard (Performance criterion B).

EN 61010-1:1993 Safety requirements for electrical equipment for measurement, control and laboratory use.

EN 61010-2-010:1995 Particular requirements for laboratory equipment for the heating of materials.

I have made all reasonable enquiries regarding the unit stated and its conformance to the following EU directives:

Low Voltage directive, 73/23/EEC and amendment 93/68/EEC, and

EMC Directive 89/336/EEC and amendments 91/263/EEC and 92/31/EEC.

To the best of my knowledge and belief these units conform to these directives.

This Declaration is controlled under an ISO 9001:2000 system certificated by BSI Quality Assurance, certificate number FM13585.

Signature

Name

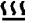
B C Coombes

Position

Quality Manager

Issue 3

23/07/2003

<p><b>Introduction</b></p> <p>Please read all the information in this booklet before using the unit.</p> <p><b>Warning</b></p> <p>HIGH TEMPERATURES ARE DANGEROUS: they can cause serious burns to operators and ignite combustible material.</p> <p>Techne have taken great care in the design of these units to protect operators from hazards, but users should pay attention to the following points:</p> <ul style="list-style-type: none"> <li>• USE CARE AND WEAR PROTECTIVE GLOVES TO PROTECT HANDS;</li> <li>• DO NOT put hot objects on or near combustible objects;</li> <li>• DO NOT operate the unit close to inflammable liquids or gases;</li> <li>• DO NOT place any liquid directly in your unit;</li> <li>• At all times USE COMMON SENSE.</li> </ul> <p><b>Operator Safety</b></p> <p>All users of Techne equipment must have available the relevant literature needed to ensure their safety.</p> <p>It is important that only suitably trained personnel operate this equipment, in accordance with the instructions contained in this manual and with general safety standards and procedures. If the equipment is used in a manner not specified by Techne the protection provided by the equipment to the user may be impaired.</p> <p>All Techne units have been designed to conform to international safety requirements and are fitted with an overtemperature cutout. On some models, the cutout is adjustable and should be set to suit the application. On all other models the cutout is preset to protect the unit.</p> <p>If a safety problem should be encountered, switch off at the mains socket and remove the plug from the supply.</p> <p><b>Installation</b></p> <ol style="list-style-type: none"> <li>1. All Techne units are supplied with a power cable. This may be integral or plug-in.</li> <li>2. Before connecting the mains supply, check the voltage against the rating plate. Connect the mains cable to a suitable plug according to the table below. <b>Note that the unit must be earthed to ensure proper electrical safety.</b></li> </ol> <table border="0"> <tr> <td><i>Connections</i></td> <td><i>220/240V</i></td> <td><i>110/120V</i></td> </tr> <tr> <td>Live</td> <td>Brown</td> <td>Black</td> </tr> <tr> <td>Neutral</td> <td>Blue</td> <td>White</td> </tr> <tr> <td>Earth</td> <td>Green/yellow</td> <td>Green</td> </tr> </table> <p>The fused plug supplied with the mains lead for use in the UK is fitted with the following value fuse to protect the cable:</p> <table border="0"> <tr> <td>230V UK</td> <td>10 AMP</td> </tr> </table> <p>The fuse in the unit protects the unit and the operator. Note that units marked 230V on the rating plate work at 220V; units marked 120V work at 110V. In both cases, however, the heating rate will degrade by approximately 8%. The rating plate is on the rear of the unit.</p> <ol style="list-style-type: none"> <li>3. Plug the mains cable into the socket on the rear of the unit.</li> <li>4. Place the unit on a suitable bench or flat workspace, or in a fume cupboard if required, ensuring that the air inlet vents on the underside are free from obstruction.</li> <li>5. Note that the following symbols may be next to the indicator lamps on the front panel of the units and have the following meanings: <ul style="list-style-type: none"> <li> : the power indicator</li> <li> : the heater indicator</li> <li> : the overtemperature indicator</li> </ul> </li> <li>6. Symbols on or near the power switch of the unit have the following meanings: <ul style="list-style-type: none"> <li>I : mains switch On</li> <li>O : mains switch Off</li> </ul> </li> </ol>	<i>Connections</i>	<i>220/240V</i>	<i>110/120V</i>	Live	Brown	Black	Neutral	Blue	White	Earth	Green/yellow	Green	230V UK	10 AMP	<p><b>After use</b></p> <p>When you have finished heating samples, remember that parts of the unit – the tubes, blocks and associated accessories – may be very hot. Take the precautions listed earlier.</p> <p><b>Guarantee</b></p> <p>The unit is guaranteed against any defect in material or workmanship for the period specified on the enclosed guarantee card. This period is from the date of purchase, and within this period all defective parts will be replaced free of charge provided that the defect is not the result of misuse, accident or negligence. Servicing under this guarantee should be obtained from the supplier.</p> <p>Notwithstanding the description and specification(s) of the units contained in the Operator's Manual, Techne hereby reserves the right to make such changes as it sees fit to the units or to any component of the units.</p> <p>This Manual has been prepared solely for the convenience of Techne customers and nothing in this Instruction Book shall be taken as a warranty, condition or representation concerning the description, merchantability, fitness for purpose or otherwise of the units or components.</p> <p><b>User maintenance</b></p> <p>NOTE THAT THIS EQUIPMENT SHOULD ONLY BE DISMANTLED BY PROPERLY TRAINED PERSONNEL. REMOVING THE SIDE, FRONT OR REAR PANELS EXPOSES POTENTIALLY LETHAL MAINS VOLTAGES. THERE ARE NO USER MAINTAINABLE PARTS WITHIN THE EQUIPMENT.</p> <p>In the unlikely event that you experience any problems with your unit which cannot easily be remedied, you should contact your supplier and return the unit if necessary. Please include any details of the fault observed and remember to return the unit in its original packing. Techne accept no responsibility for damage to units which are not properly packed for shipping: if in doubt, contact your supplier. See the Decontamination Certificate supplied with your unit.</p> <ol style="list-style-type: none"> <li>1. <b>Cleaning</b></li> </ol> <p>Before cleaning your unit ALWAYS disconnect from the power supply and allow to cool below 50° C.</p> <p>Your unit can be cleaned by wiping with a damp soapy cloth. Care should be exercised to prevent water from running inside the unit. Do not use abrasive cleaners.</p> <ol style="list-style-type: none"> <li>2. <b>Fuses</b></li> </ol> <p>Your unit is protected by one or two fuses. These should only be changed by suitably qualified personnel. If the fuses blow persistently, a serious fault is indicated and you may need to return the unit to your supplier for repair.</p> <p><b>Contact Information</b></p> <p>For technical, sales or servicing information, contact your local Techne dealer or,</p> <p>Techne, Duxford, Cambridge, CB2 4PZ, United Kingdom.</p> <p>Telephone: +44(0)1223 832401  Fax: +44(0)1223 836838  Service: +44(0)1223 836950 Out of office hours  e-mail: sales@techne.com  Web site: www.techne.com</p> <p>or,</p> <p>Techne Inc, 3 Terri Lane, Suite 10, Burlington, New Jersey 08016, USA.</p> <p>Telephone: 609-589-2560  Toll free: 800-225-9243 ext 306  Fax: 609-589-2571  e-mail: labproducts@techneusa.com  Web site: www.techneusa.com</p>
<i>Connections</i>	<i>220/240V</i>	<i>110/120V</i>													
Live	Brown	Black													
Neutral	Blue	White													
Earth	Green/yellow	Green													
230V UK	10 AMP														

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## THE CALIBRATOR

The Tecal Calibrator is designed to provide safe and convenient calibration of a wide range of thermal sensors. It features fast heat up times, with accuracy and repeatability.

### Applications

**Before using the Tecal, make sure you have read this manual carefully. If you have any queries, contact your supplier.**

The Tecal unit can calibrate temperature probes without the need to return them to a specialist laboratory. To ensure accuracy the unit must be operated in an environment with a stable ambient temperature.

The thermal sensors are placed in a well in the temperature controlled block. A number of inserts with different dimension wells are available to match standard probe sizes: these are detailed towards the back of this manual.

The required temperature is set on the calibrator and the operation of the probe is checked.

A UKAS certificate can be provided, contact your dealer for details.

### Warning

**HIGH TEMPERATURES ARE DANGEROUS:** they can cause serious burns to operators and ignite combustible material.

Techne have taken great care in the design of these units to protect operators from hazards, but operators should pay attention to the following points:

- USE CARE AND WEAR PROTECTIVE GLOVES TO PROTECT HANDS;
- DO NOT put hot objects on or near combustible objects;
- DO NOT operate the unit close to inflammable liquids or gases;
- DO NOT place any liquid directly in your Tecal unit;
- DO NOT place the unit into its carrying case if the block temperature exceeds 50°C;
- If you are using a cooling probe, make sure that water is flowing through the cooling probe BEFORE inserting the probe into the block;
- At all times USE COMMON SENSE.

## Specification

**The figures stated in the specification are as measured 0-50 mm from the base of the well of the insert, using an SPRT, at the time of calibration.**

Measuring zone	0 to 50 mm	0" to 2"
	from base of well	from base of well
Depth of well in insert	114.3 mm	4.5"
Temperature range	20°C above ambient to 425°C	36°F above ambient to 797°F
Temperature accuracy in measuring zone	±0.3°C	±0.54°F
Temperature uniformity in measuring zone	±0.2°C at 300°C	0.36°F at 542°C
Temperature stability after 10 min	±0.05°C	±0.09°F
Display resolution	0.1°C	0.1°F
Heat up rate	15 minutes	20°C to 400°C 68°F to 752°F
Cool down rate	25 minutes	400°C to 100°C 752°F to 182°F
Fan cooling	Automatic	
Comms Port 9way D type	Full bi-directional RS232	

### Dimensions

Height	272mm	10.7"
Width	185mm	7.28"
Depth	439mm	17.28"

### Wells in the block:

Central Well for insert, see the list of inserts at the back of this book

Depth	130mm	5.1"
Diameter	38.2mm	1.5"

Reference well\*

Depth	101.6mm	4"
Diameter	6.6mm	0.26"

## Working conditions

The Tecal units are designed to work safely under the following conditions:

Ambient temperature range 5°C to 40°C

Humidity Up to 95% relative humidity, non-condensing

**Note:** The control specifications quoted are for an ambient temperature range of 10°C to 30°C. The specification may deteriorate outside this range but the unit will still work safely.

Radio frequency interference tested and passed to EN50081-1.


Immunity Tested and passed to EN50082-1

## Packing

When you receive your unit, make sure you keep the original packing in case you ever need to return it for service or repair. When returning a unit, remember to remove the insert from the temperature controlled block.


The unit must be transported in the original packing to avoid damage. The packing comprises; the unit in the soft carrying case packed with foam into an outer cardboard box. Tecal accepts no responsibility for damage incurred unless the unit is correctly packed and transported in this way.

When you first switch on, the display will show the edition of the software which your unit uses. For example software issue "1" would be shown as follows

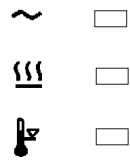
It will display this for 1 second, then the actual temperature of the block will be indicated.


The front panel controls consist of three buttons for controlling the display, a five digit LED display and three indicators.

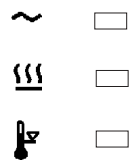The **SET** temperature button displays the set temperature when pressed.

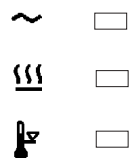
When the **SET** temperature button is held down and the **UP ARROW** button is pressed, the set temperature is increased.

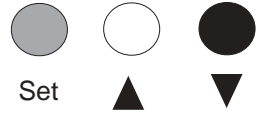When the **SET** temperature button is held down and the **DOWN ARROW** button is pressed, the set temperature is decreased.

**Speed of Change of Set Temperature**

Each press of the **UP ARROW** or **DOWN ARROW** buttons will increase or decrease the set temperature by 0.1°C. If the buttons are held down the temperature change will accelerate to 5° per second

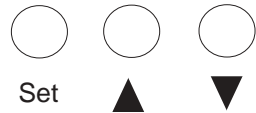The top indicator shows that there is power to the unit

The next indicator shows when there is power to the block.

When the set temperature is higher than the block temperature this indicator will light. If the light is on continuously the heater is getting constant power. The only exception is described under Over-Temperature Indicator. As the temperature approaches the set temperature the indicator will flash. When set temperature is reached the indicator will stay on for shorter periods. If the block temperature is above the set temperature then the indicator will be off, as the heater is not getting any power.

If the unit should, for any reason, exceed the temperature set for the over-temperature cutout (see page 15), the over-temperature indicator will light. The heater will have been switched off and the unit will begin to cool even if the heater light is on (the light staying on or not depends on which circuit has sensed an over-temperature).

If the unit there should, for any reason, be a sensor fault, the bottom indicator will light. The power to the block will have been switched off and the unit will begin to return to ambient even if the power light is on (the light staying on or not depends on which circuit has sensed a fault).

### Setting the units for temperature

The unit will normally be set to display in °C, if you require it to work and display in °F press both the **UP ARROW** and **DOWN ARROW** buttons at the same time. To convert back, press both buttons again.

### Setting the Operating Temperature

1. To display the set temperature on the digital display, press and hold the **SET** Temperature button.  
To adjust the set temperature, press the **SET** Temperature button and hold it while pressing the **UP** or **DOWN** buttons. When the **SET** Temperature button is released, the measured temperature is displayed (in degrees Celsius).
2. The heater (and heater indicator) comes on if the set temperature is higher than the current block temperature. The fan will automatically come on if the set temperature is more than 15°C below the actual temperature of the block.
3. When the measured temperature approaches the set temperature, the heater indicator will begin to flash. As the measured temperature stabilises the indicator will stay on for shorter periods.
4. There will be a time-lag between the heater block and the insert achieving the set temperature due to thermal contact between them. The display may reach temperature in 2 or 3 minutes while the insert may take 5 minutes or more to reach temperature.

## Cooling Probe

Operating temperatures below ambient can be achieved by using the cooling probe with a chilled water circulator. The cooling probe is supplied as an accessory and can be used in an insert with a 3/8" or 10mm hole (see Accessories). Instructions for its use are included with the probe.

The block should not be operated at temperatures below dewpoint as excessive build up of moisture will occur which may cause damage to the interior of the unit.

## After Use

1. When you have finished heating samples, remember that parts of the unit – inserts and associated accessories – may be very hot. Take the precautions listed earlier. We recommend that the inserts should be allowed to cool to 70°C before being removed from the Tecal unit. They will still have to be handled with care
2. If you need to remove an insert while it is hot, fit the extractor tool into the locating holes and lift the insert out carefully. Never leave the extractor tool in the insert while it is being used in the Tecal unit.


## RS232 SERIAL INTERFACE

The calibrator may send data logging information to an IBM PC or compatible computer by connecting the unit and the PC via an RS232 cable, and installing the "Calsoft" software. Contact your supplier for details.

**Ensure that there is a connection made before the calibrator is switched on.**

The RS232 cable must be fitted to both the unit and the PC before either unit is powered up, otherwise, data integrity cannot be guaranteed. Once the cable is fitted, it does not matter which unit is powered up first.

RS232 9-way PC serial port

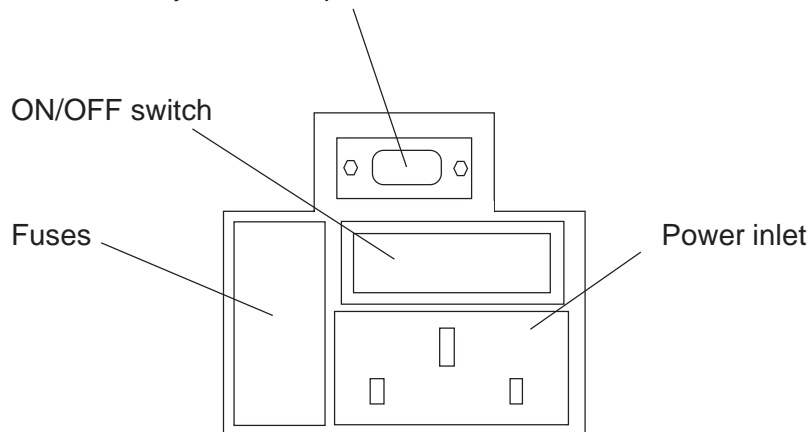The following tables indicate the cable specifications for a 9-way PC serial port:

<b>Calibrator</b>		<b>PC</b>
<i>9-way female D type</i>		<i>9-way female D type</i>
pin	signal	pin
Case	F.GND	Case
3	TxD	3
2	RxD	2
7	RTS	7
8	CTS	8
6	DSR	6
1	CD	1
4	DTR	4
5	S.GND	5

### Sending data to the PC

The procedure for sending data to the PC is described in the instructions supplied with the software. The PC must be properly connected by the RS232 cable and running Calsoft to accept or display data.

## TECHNICAL INFORMATION

Brief fault finding notes and lists of replacement parts, accessories and inserts for the Tecal are given in this section.

NOTE THAT THIS EQUIPMENT SHOULD ONLY BE DISMANTLED BY PROPERLY TRAINED PERSONNEL. **REMOVING THE OUTER COVERS OR BASE EXPOSES POTENTIALLY LETHAL MAINS VOLTAGES.**

THERE ARE NO OPERATOR SERVICEABLE PARTS WITHIN THE EQUIPMENT.

### General advice

In the unlikely event that you experience any problems with your Tecal which cannot be easily remedied, you should contact your supplier and return the unit if necessary. Please include any details of the fault observed and remember to return the unit in its original packing. The insert must be removed from the unit and packed separately within the case. Techne accept no responsibility for damage to units which are not properly packed for shipping: if in doubt, contact your supplier.

### General fault finding

#### Over temperature cutout

Your Tecal unit is fitted with two independent circuits to protect it from overheating. The unit constantly checks that the block temperature does not exceed its maximum. If for some reason this temperature is exceeded, all power to the block is cut.

Allow the unit to cool to a safe temperature (less than 100°C) before switching off and check for any obvious causes of overheating before switching it back on. Switching off the mains power resets the over-temperature cutout.

Repeated cutouts indicate a serious fault and you should return the unit to your supplier for repair.

#### Cleaning your Tecal unit

Before cleaning your unit, disconnect from the power supply and allow to cool to ambient temperature.

You can clean the case of the Tecal with a cloth dipped in water or ethanol (methanol can also be used). No part of the case or cover should be immersed in the solvents.

Do not use acetone or abrasive cleaners.

Before using any cleaning or decontamination method except those recommended in this manual, the responsible body should check with Techne that the proposed method will not damage the equipment.

#### HELP display

If HELP is displayed there is a problem with the internal sensing device. The calibrator should be returned to your supplier for repair.

**Fuses**

If neither the power light nor display (on the front panel) is lit, one of the two fuses may have blown. Make sure there is no external cause (such as a faulty plug or lead).

Fuses should only be changed by suitably qualified personnel. If the fuses blow persistently, a serious fault is indicated and you should return the calibrator to your supplier for repair.

Never fit a fuse rated higher than the value indicated on the unit, serious damage or personal injury may result.

**Replacement parts**

Each unit is supplied with an extraction tool, a soft carrying case and a mains cable. The following parts may be purchased if replacements or alternatives are required.

6101058	Insert extractor
FCAB10UK	Mains cable, 230V UK units
FCABLEEU	Mains cable, 230V units
FCABLEUS	Mains cable, 120/100V units
6103711	Soft carrying case

**Accessories**

<i>Part N°</i>	<i>Description</i>
FDB00CP	Cooling Probe
FCAL232	CALSOFT PC software kit

## Inserts

Inserts are made from aluminium or aluminium-bronze and must be ordered separately from the calibrator. Each insert is stamped for identification.

A blank insert is available which may be drilled to suit your application.

Immersion depth 114.3 mm or 4.5"

*NUMBER*      *PROBE DIAMETER*

### **METRIC**

FINSALA	5 off 6mm
FINSALB	1 off each 10mm, 8mm, 6mm, 4.5mm, 3mm
FINSALC	2 off 6mm, 2 off 10mm
FINSALD	2 off 6mm, 2 off 12mm
FINSALE	1 off 6mm
FINSALZ	1 off 20mm

### **IMPERIAL**

FINSALF	5 off 1/4"
FINSALG	1 off each 3/8", 5/16", 1/4", 3/16", 1/8"
FINSALH	2 off 1/4", 2 off 3/8"
FINSALI	2 off 1/4", 2 off 1/2"
FINSALJ	1 off 1/4"
FINSALL	1 off 9/16"
FINSALM	1 off 5/8"
FINSALN	1 off 3/4"
FINSALO	1 off 11/16"

### **BLANK**

FINSALK	BLANK
---------	-------